

Mini-spectrometers



[TG series]

C11713CA

C11714CB

**For Raman spectrophotometry,
high resolution type (spectral resolution: 0.3 nm)**

The C11713CA and C11714CB are polychromators integrated with optical elements, an image sensor and a driver circuit. Light to be measured is guided into the entrance port of TG series through an optical fiber and the spectrum measured with the built-in image sensor is output from the USB port to a PC for data acquisition. The C11713CA has sensitivity in a wavelength range of 500 to 600 nm, while the C11714CB covers a range of 790 to 920 nm. Both types offer a spectral resolution of 0.3 nm. Their products come supplied with free evaluation software that allows setting measurement conditions, acquiring and saving data, and displaying graphs. Original measurement software can be designed on an end-user's side as DLL's function specification is disclosed.

Features

- High accuracy optical characteristics: spectral resolution 0.3 nm
- Easy to install into equipment due to compact design
- High throughput due to transmission grating made of quartz
- Wavelength conversion factor is recorded in internal memory *1

*1: A conversion factor for converting the image sensor pixel number into a wavelength is recorded in the module. A calculation factor for converting the A/D converted count into the input light intensity is not provided.

Applications

- Raman spectrophotometry


Selection guide (Typ.)

Parameter	C11713CA	C11714CB	Unit
Spectral response range	500 to 600	790 to 920	nm
Spectral resolution	0.3		nm
Built-in image sensor	Back-thinned CCD image sensor	IR-enhanced back-thinned CCD image sensor	-

High sensitivity in near infrared region (C11714CB)

Using a CCD image sensor with enhanced near infrared sensitivity, the C11714CB offers higher sensitivity than that of the previous product.

Spectral response

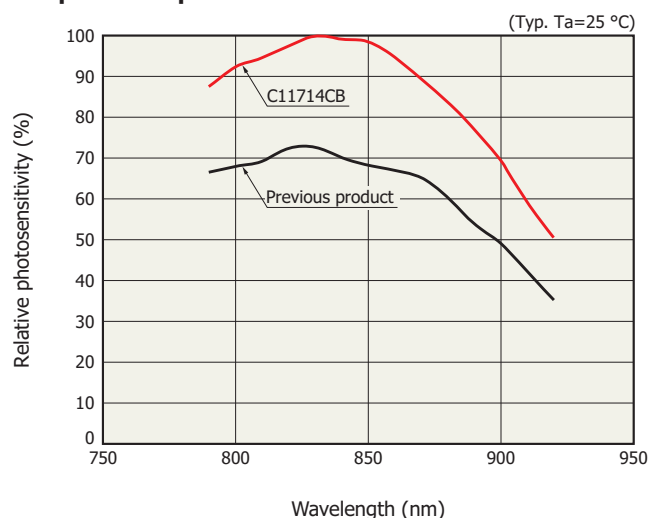

KACC80371EA

Structure

Parameter	C11713CA	C11714CB	Unit
Dimensions (W × D × H)	120 × 70 × 60		mm
Weight	592		g
Image sensor	Back-thinned type CCD image sensor (S10420-1106-01)	IR-enhanced back-thinned CCD image sensor (S11510-1006)	-
Number of pixels	2048	1024	pixels
Slit*2 (H × V)	10 × 1000		μm
NA*3	0.11		-
Connector for optical fiber	SMA905D		-
A/D conversion	16		bit
Interface	USB 1.1		
External power supply (AC adapter)	AC 100V-240V, 50/60Hz		°C

*2: Entrance slit aperture size

*3: Numeric aperture (solid angle)

Absolute maximum ratings

Parameter	Value	Unit
External power supply voltage	5.25	V
Operating temperature*4	+5 to +40	°C
Storage temperature*4	-20 to +70	°C

*4: No dew condensation

When there is a temperature difference between a product and the surrounding area in high humidity environment, dew condensation may occur on the product surface. Dew condensation on the product may cause deterioration in characteristics and reliability.

Note: Exceeding the absolute maximum ratings even momentarily may cause a drop in product quality. Always be sure to use the product within the absolute maximum ratings.

Electrical characteristics (Ta=25 °C, unless otherwise noted)

Parameter	Min.	Typ.	Max.	Unit
Integration time	10 to 10000			ms
Consumption current of USB bus power*5	-	-	150	mA
External power supply voltage	-	5	-	V
Consumption current of external power supply*5	-	-	0.8	A

*5: When operated with the supplied evaluation software (at default settings, dark state, excluding start-up).

Electrical and optical characteristics (Ta=25 °C, unless otherwise noted)



Parameter	C11713CA			C11714CB			Unit
	Min.	Typ.	Max.	Min.	Typ.	Max.	
Spectral response range	500 to 600			790 to 920			nm
Spectral resolution (FWHM)*6	-	0.3	0.5	-	0.3	0.5	nm
Wavelength reproducibility*7	-0.1 to +0.1			-0.1 to +0.1			nm
Wavelength temperature dependence	-0.04 to +0.04			-0.04 to +0.04			nm/°C
Spectral stray light*6 *8	-	-	-30	-	-	-30	dB

*6: Depends on the slit opening. Values were measured with the slit listed in the table "Structure".



*7: Measured under constant light input conditions

*8: When monochromatic light of the following wavelengths is input, spectral stray light is defined as the ratio of the count measured at the input wavelength, to the count measured in a region of the input wavelength ±10 nm
C11713CA: 550 nm, C11714CB: 860 nm

▣ Spectral response (typical example)

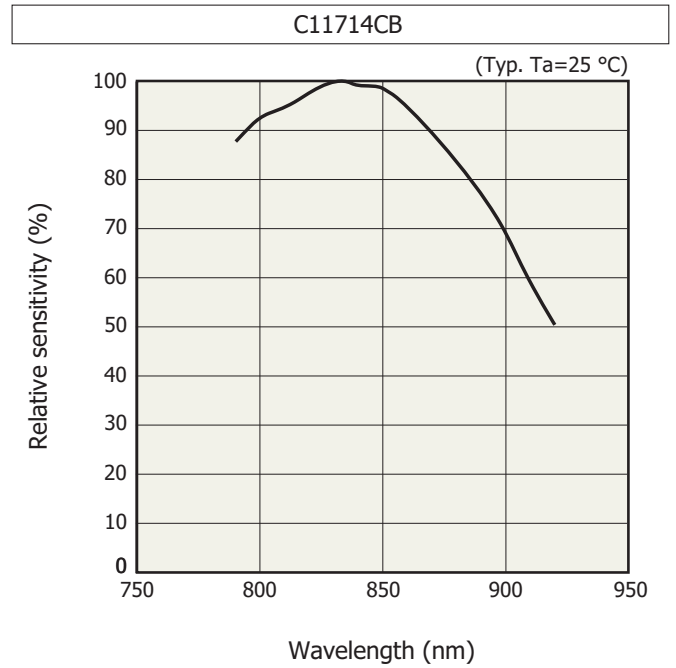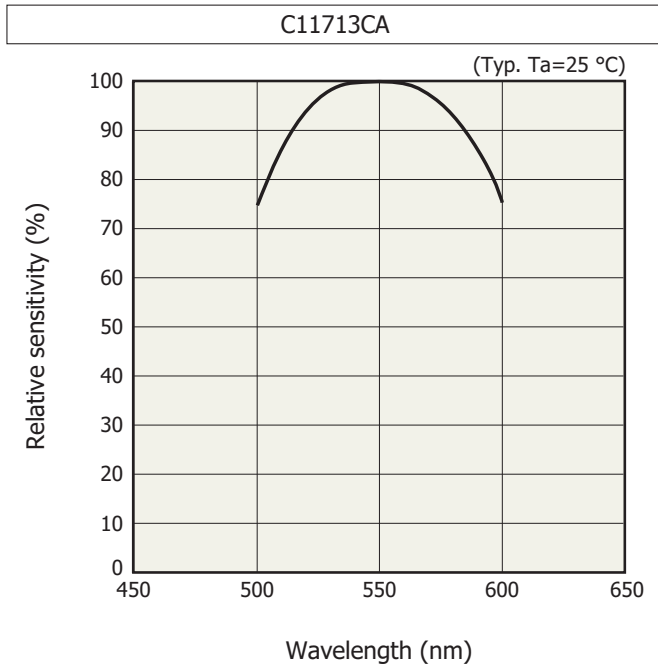


▣ Spectral resolution vs. wavelength

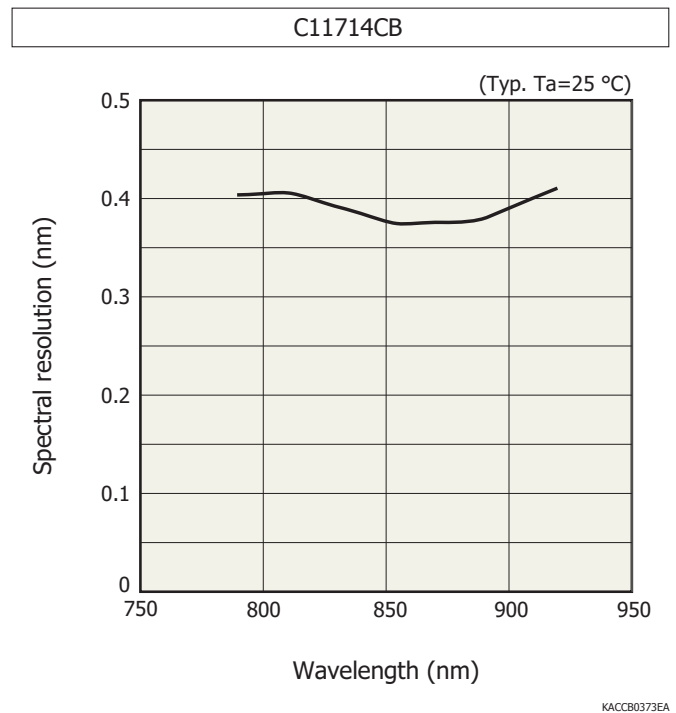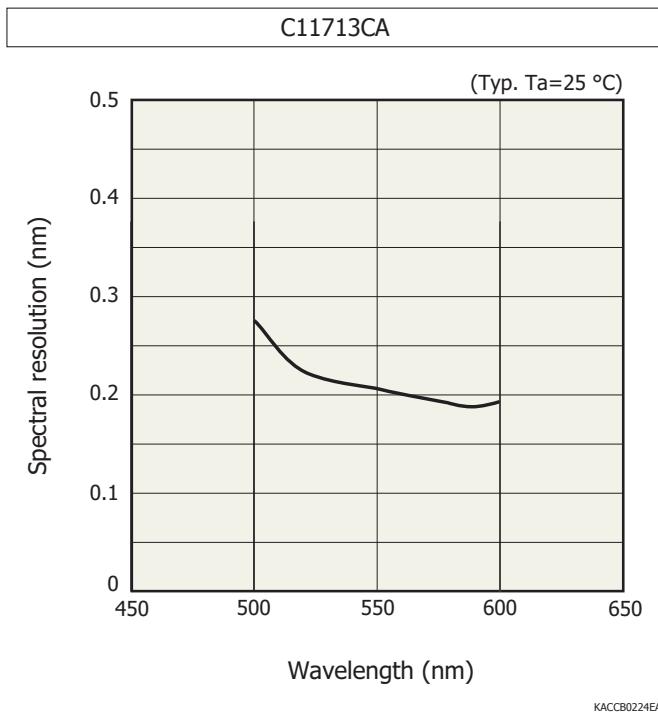


Linearity (typical example)

C11713CA

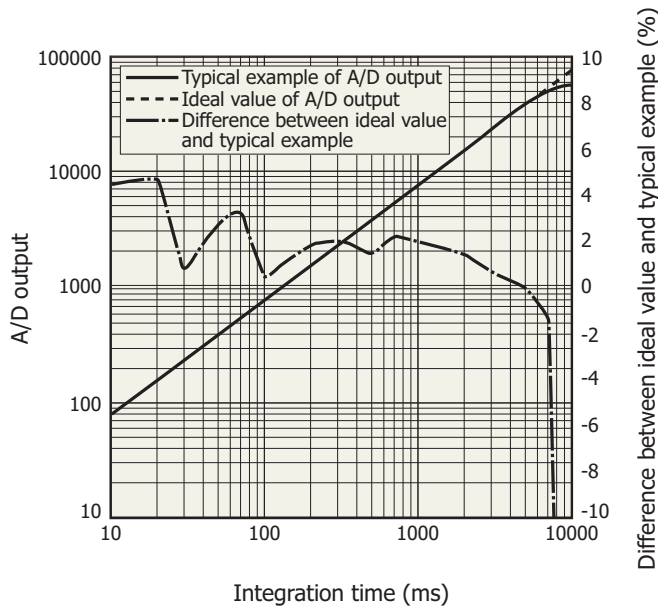

KACCB0235EA

C11714CB

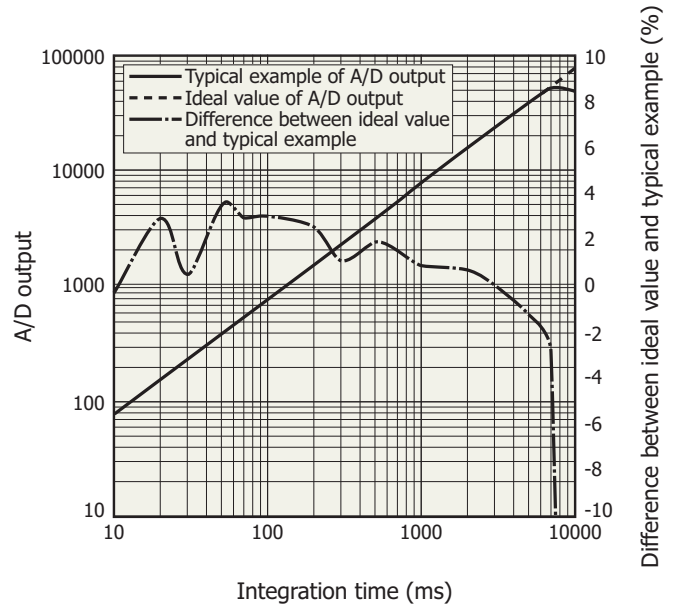


KACCB0374EA

A/D output is the output with dark output is subtracted when light is input. The difference between the ideal value and typical example contains a measurement error. The smaller the A/D output, the larger the measurement error.


Dark output vs. ambient temperature (typical example)

C11713CA

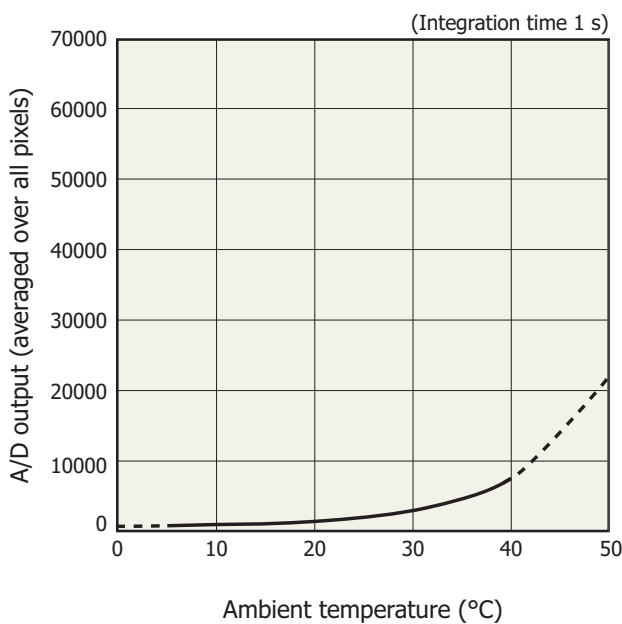

KACCB0237EA

C11714CB

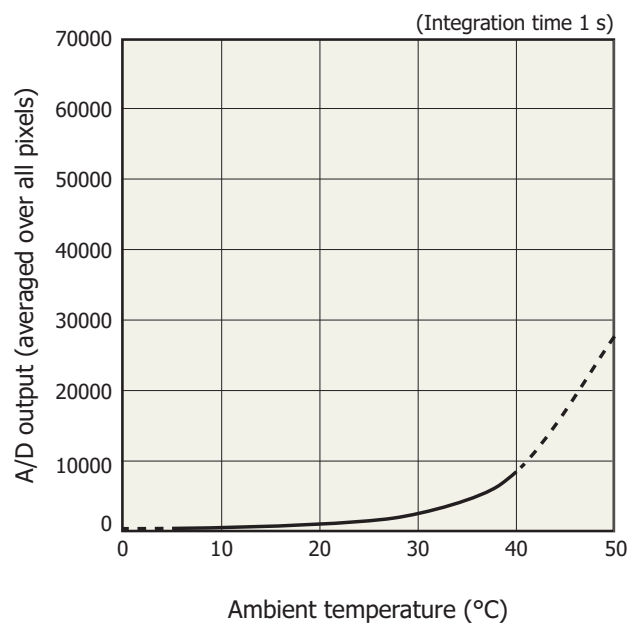


KACCB0375EA

A/D output is the sum of the sensor and circuit offset outputs and the sensor dark output.


Connection example

Light to be measured is guided into the entrance port of TG series through an optical fiber and the spectrum measured with the built-in image sensor is output through the USB port to a PC for data acquisition. There are no moving parts inside the unit so stable measurements are obtained at all times. An optical fiber that guides light input from external sources allows a flexible measurement setup.

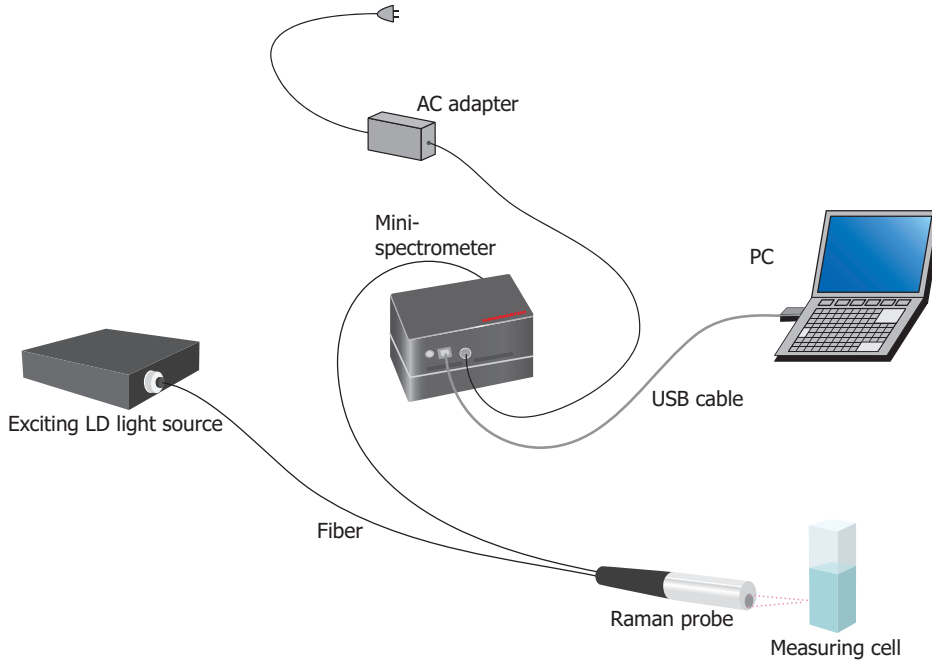


Optical component layout

TG series mini-spectrometers use a transmission holographic grating made of quartz and precision optical components arranged on a rugged optical base, making it possible to deliver high throughput and highly accurate optical characteristics.

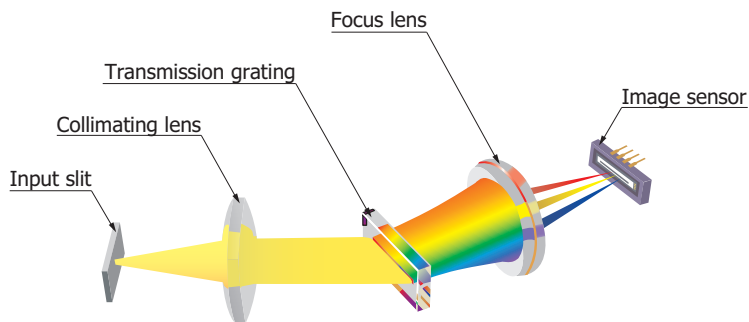

Measurement example

Raman spectrophotometry using nophthalene (C11714CB)

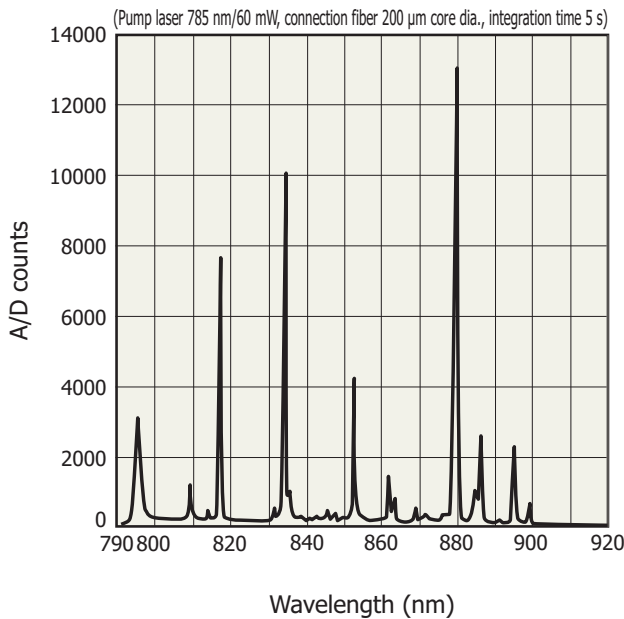

KACCB0229EA

Evaluation software package (supplied with unit)

Installing the evaluation software package (Spec Evaluation.exe)*8 into your PC allows running the following basic tasks:

- Measurement data acquisition and save
- Measurement condition setup
- Module information acquisition
(wavelength conversion factor, polychromator type, etc.)
- Graphic display
- Arithmetic operation
 - Pixel number to wavelength conversion
 - Comparison calculation with reference data
(transmittance, reflectance)
 - Dark subtraction
 - Gaussian approximation
(peak position and count, FWHM)

Note:

- Two or more mini-spectrometers can be connected and used with one PC simultaneously.
- The external trigger input function does not work with the evaluation software. If using an external trigger input or designing original application software, the user software must be configured to support that function.

*8: Compatible OS: Microsoft® Windows® 7 Professional SP1 (32-bit, 64-bit)
 Microsoft Windows 8 Professional (32-bit, 64-bit)
 Microsoft Windows 10 Professional (32-bit, 64-bit)


DLL for controlling hardware is also provided.

You can develop your own measurement programs by using a following software development environment.


Microsoft Visual Studio® 2008 (SP1) Visual C++®

Microsoft Visual Studio 2008 (SP1) Visual Basic®

Note: Microsoft, Windows, Visual Studio, Visual C++ and Visual Basic are either registered trademarks or trademarks of Microsoft Corporation in the United States and other countries.

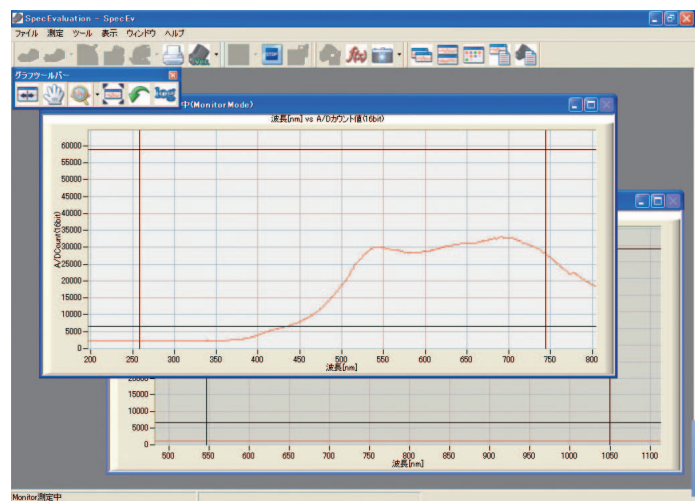

Dimensional outline (unit: mm, tolerance unless otherwise noted: ±0.5)

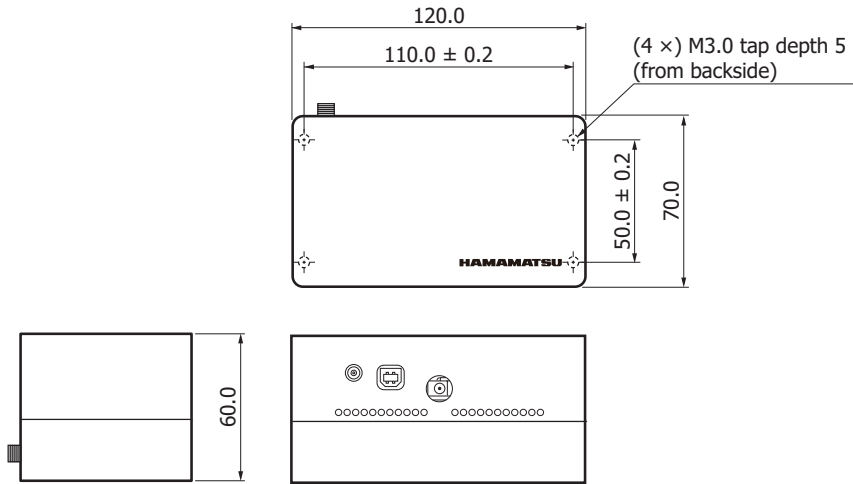

Weight: 592 g

KACCA0281EB

Accessories


- USB cable
- Dedicated software (evaluation software, sample software, DLL)
- AC adapter (for power supply)

Options (sold separately)

- Coaxial cable for external trigger input A10670
- Optical fibers for light input

Type no.	Product name	Core diameter (μm)	Specification
A9763-01	Fiber for visible/near infrared range	600	NA=0.22, length 1.5 m, connectorized SMA905D at both ends

Dimensional outline (A10670, unit: mm)

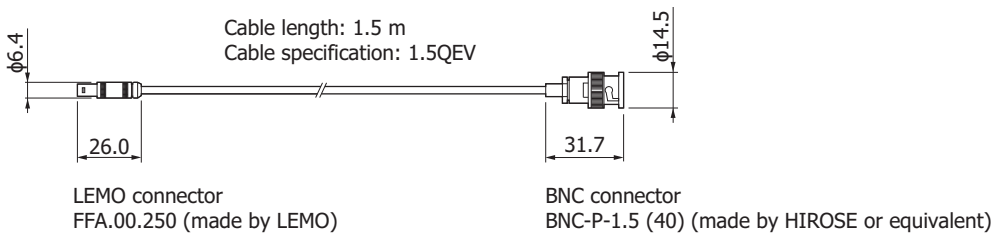

KACCA0220EB

Mini-spectrometer/micro-spectrometer lineup

Type no.	Type	Spectral response range (nm)													Spectral resolution max. (nm)	Image sensor		
		200	400	600	800	1000	1200	1400	1600	1800	2000	2200	2400	2600				
C10082CA	TM-UV/VIS-CCD High sensitivity																6	Back-thinned CCD image sensor
C10082CAH	TM-UV/VIS-CCD High resolution		200 to 800														1*	
C10082MD	TM-UV/VIS-MOS Wide dynamic range																6	CMOS linear image sensor
C10083CA	TM-VIS/NIR-CCD High sensitivity																8 (λ=320 to 900 nm)	Back-thinned CCD image sensor
C10083CAH	TM-VIS/NIR-CCD High resolution																1* (λ=320 to 900 nm)	
C10083MD	TM-VIS/NIR-MOS Wide dynamic range																8	CMOS linear image sensor
C11697MB	TM-VIS/NIR-MOS-II Trigger-compatible																8	High-sensitivity CMOS linear image sensor
C9404CA	TG-UV-CCD High sensitivity		200 to 400														3	Back-thinned CCD image sensor
C9404CAH	TG-UV-CCD High resolution																1*	
C9405CB	TG-SWNIR-CCD-II IR-enhanced																5 (λ=550 to 900 nm)	IR-enhanced back-thinned CCD image sensor
C11713CA	TG-RAMAN-I High resolution																0.3*	Back-thinned CCD image sensor
C11714CB	TG-RAMAN-II High resolution																0.3*	IR-enhanced back-thinned CCD image sensor
C11482GA	TG2-NIR Non-cooled type																7	InGaAs linear image sensor
C9913GC	TG-cooled NIR-I Low noise (cooled type)																7	
C9914GB	TG-cooled NIR-II Low noise (cooled type)																8	
C11118GA	TG-cooled NIR-III Low noise (cooled type)																20	
C13053MA	TF-SWIR-MOS-II Compact, thin case																3.5	High-sensitivity CMOS linear image sensor
C13054MA	TF-RAMAN Compact, thin case																0.4*	
C14214MA	TF-RAMAN Compact, thin case																0.6	
C13555MA	TF-VIS-MOS-II Compact, thin case																3	
C11007MA	RC-VIS-MOS Spectrometer module																9	CMOS linear image sensor
C11008MA	RC-SWNIR-MOS Spectrometer module																8	IR-enhanced CMOS linear image sensor

* Typ.

For installation into mobile measuring equipment

Type no.	Type	Spectral response range (nm)													Spectral resolution max. (nm)	Image sensor		
		200	400	600	800	1000	1200	1400	1600	1800	2000	2200	2400	2600				
C11009MA	RC-VIS-MOS Spectrometer head																9	CMOS linear image sensor
C11010MA	RC-SWNIR-MOS Spectrometer head																8	IR-enhanced CMOS linear image sensor

For installation into mobile measuring equipment (ultra-compact)

Type no.	Type	Spectral response range (nm)													Spectral resolution max. (nm)	Image sensor		
		200	400	600	800	1000	1200	1400	1600	1800	2000	2200	2400	2600				
C11708MA	MS-SWNIR-MOS Spectrometer head																20	CMOS linear image sensor
C12666MA	Spectrometer head																15	CMOS linear image sensor
C12880MA	Spectrometer head																15	High-sensitivity CMOS linear image sensor

Related information

www.hamamatsu.com/sp/ssd/doc_en.html

■ Precautions

- Disclaimer
- Mini-spectrometers

■ Technical information

- Mini-spectrometers

Information described in this material is current as of July 2018.

Product specifications are subject to change without prior notice due to improvements or other reasons. This document has been carefully prepared and the information contained is believed to be accurate. In rare cases, however, there may be inaccuracies such as text errors. Before using these products, always contact us for the delivery specification sheet to check the latest specifications.

The product warranty is valid for one year after delivery and is limited to product repair or replacement for defects discovered and reported to us within that one year period. However, even if within the warranty period we accept absolutely no liability for any loss caused by natural disasters or improper product use. Copying or reprinting the contents described in this material in whole or in part is prohibited without our prior permission.

HAMAMATSU

www.hamamatsu.com

HAMAMATSU PHOTONICS K.K., Solid State Division

1126-1 Ichino-cho, Higashi-ku, Hamamatsu City, 435-8558 Japan, Telephone: (81) 53-434-3311, Fax: (81) 53-434-5184

U.S.A.: Hamamatsu Corporation: 360 Foothill Road, Bridgewater, N.J. 08807, U.S.A., Telephone: (1) 908-231-0960, Fax: (1) 908-231-1218, E-mail: usa@hamamatsu.com

Germany: Hamamatsu Photonics Deutschland GmbH: Arzbergerstr. 10, D-82211 Herrsching am Ammersee, Germany, Telephone: (49) 8152-375-0, Fax: (49) 8152-265-8, E-mail: info@hamamatsu.de

France: Hamamatsu Photonics France S.A.R.L.: 19, Rue du Saule Trapu, Parc du Moulin de Massy, 91882 Massy Cedex, France, Telephone: 33-(1) 69 53 71 00, Fax: 33-(1) 69 53 71 10, E-mail: infos@hamamatsu.fr

United Kingdom: Hamamatsu Photonics UK Limited: 2 Howard Court, 10 Tewin Road, Welwyn Garden City, Hertfordshire AL7 1BW, United Kingdom, Telephone: (44) 1707-294888, Fax: (44) 1707-325777, E-mail: info@hamamatsu.co.uk

North Europe: Hamamatsu Photonics Norden AB: Torshamnsgatan 35 16440 Kista, Sweden, Telephone: (46)8-509 031 00, Fax: (46)8-509 031 01, E-mail: info@hamamatsu.se

Italy: Hamamatsu Photonics Italia S.r.l.: Strada della Moia, 1 int. 6, 20020 Arese (Milano), Italy, Telephone: (39)02-93 58 17 33, Fax: (39)02-93 58 17 41, E-mail: info@hamamatsu.it

China: Hamamatsu Photonics (China) Co., Ltd.: B1201, Jiaming Center, No.27 Dongsanhuan Beilu, Chaoyang District, Beijing 100020, China, Telephone: (86) 10-6586-6006, Fax: (86) 10-6586-2866, E-mail: hpc@hamamatsu.com.cn

Taiwan: Hamamatsu Photonics Taiwan Co., Ltd.: 8F-3, No. 158, Section2, Gongdao 5th Road, East District, Hsinchu, 300, Taiwan R.O.C. Telephone: (886)03-659-0080, Fax: (886)03-659-0081, E-mail: info@hamamatsu.com.tw