



Near Infrared Oxygenation Monitor

NIRO-200NX



HAMAMATSU

A highly functional tissue oxygenation monitor to meet a variety of needs

- Continuous, non-invasive measurement of tissue oxygenation using low light.
- Easy operation: simply attach the reusable probes and press the start button.
- Useful in a wide variety of applications, from management of brain oxygenation in operating rooms to clinical studies related to brain function, brain metabolism, muscle function, etc.

Advantages of measuring with the NIRO-200NX

- Measurement using both MBL (ΔO_2Hb , ΔHHb , ΔcHb) and SRS (TOI and THI) methods, running independently of each other
- Combination of MBL and SRS can provide local haemodynamic information such as changes in oxygen consumption or vasomotor activities, and provides increased data reliability.
- Attachments are reusable.

Applications:

Monitoring brain oxygenation during heart surgery and other procedures

Patient monitoring in emergency rooms

Patient monitoring in ICUs, CCUs, and NICUs

Clinical studies related to oxygenation levels and blood metabolism in the brain

Studies of muscle tissue oxygenation for sports medicine, rehabilitation and similar applications

Can be used on emergency carts

Can be used on rescue helicopters



Highly improved features for clinical use


Improved features:

Compact and lightweight

Half the size of the previous NIRO-200, providing the same great performance while making it much easier to use with an emergency cart or in the operating theatre.

External monitor output

View data on an external monitor at the same time, useful for anaesthetists and perfusionists

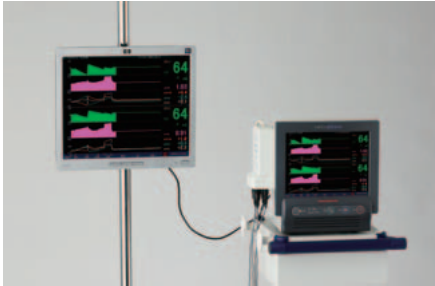

Past data retrievable

Data stored in the NIRO-200NX can be retrieved for purposes such as patient counselling or research.

Touch panel display

The machine can be operated by touching the screen.



Input of event marks from external devices

Ability to input event marks externally, useful for analyzing the results of pump control and for research purposes

Equipped with battery

Useful in case of power failure or during patient transfer after an operation


Data storage to USB memory devices

Data may be saved to a USB memory device, useful for analyzing the results of pump control.



Basic features of the NIRO-200NX

The NIRO-200NX is a tissue oxygenation monitor that uses near infrared spectroscopy. The unit uses safe, low light to measure the Tissue Oxygenation Index (TOI), showing the oxygen saturation level, the Normalised Tissue Haemoglobin Index (nTHI), showing the percentage change in the amount of initial haemoglobin, as well as changes in concentration of oxygenated haemoglobin (ΔO_2Hb), deoxygenated haemoglobin (ΔHHb), and total haemoglobin (ΔcHb), all in real time.

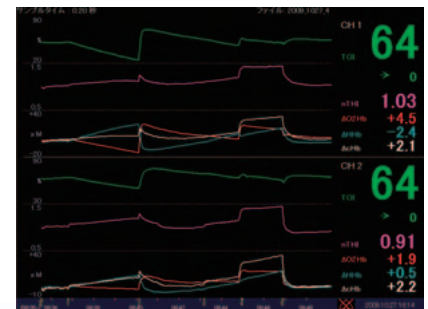

- Measure data simultaneously in two locations (two channels installed)
- Measure concentration changes in oxygenated, deoxygenated and total haemoglobin
- Easy operation by simply pressing the start button
- Transfer numerical value data to a personal computer
- Safety and high reliability



NIRO-200NX

Specification:


Measurement items	Tissue oxygenation index, TOI (%) Normalized tissue hemoglobin index, nTHI (absolute value in arbitrary unit) Oxygenated hemoglobin change, ΔO_2Hb ($\mu\text{mol/L}$) Deoxygenated hemoglobin change, ΔHHb ($\mu\text{mol/L}$) Total hemoglobin change, ΔcHb ($\mu\text{mol/L}$)
Sample interval	0.05 s, 0.2 s, 0.5 s, 1 s, 2 s, 5 s, 10 s, 20 s, 30 s
Light source	LED (735 nm, 810 nm, 850 nm: nominal values)
Output power	Less than 2 mW
Light detector	Photodiode
Measurement method	SRS method (Spatially Resolved Spectroscopy) and MBL method (Modified Beer-Lambert)
Data memory	Internal backup memory for more than 10 data files. (Maximum number of data samples for each file is 50 000.)
Saving data to a USB device	Measured data can be saved to an external USB memory device using the USB connector.
Output signal	Digital output (RS-232C) / Analog output / Philips Vuelink format
Measurement probes	Approx. 2.5 m cable length (for emission and detection probes)
Battery	Operating time of about 30 min. (fully charged)
External Event Input	TTL, "L" with a duration longer than 100 ms
Display Unit	Size: approx. 264 mm(W) \times 279 mm(H) \times 190 mm(D) Weight: approx. 6.0 kg
AMP Unit	Size: 91 mm(W) \times 55 mm(H) \times 156 mm(D) Weight: approx. 0.7 kg Cable length: approx. 4 m
Power consumption	Less than 90 V·A

Standard configuration:

· Display Unit (DU).....	1
· AMP Unit (AU).....	1
· Emission probe (A10959).....	2
· Detection probe (A10962).....	2
· Probe holder S type (A10963).....	2
· Probe holder L type (A10965).....	2
· Double-sided adhesive tape for S type (A10967).....	1 bag
· Double-sided adhesive tape for L type (A10968).....	1 bag
· Cable clip (A11047).....	1 bag
· Power supply cord.....	1
· Manual.....	1

Options:

- Data download software, Windows compatible (U10898)

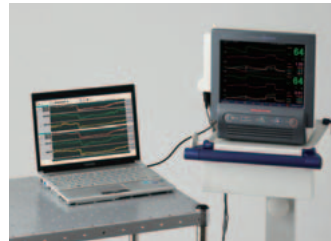

The data can be displayed and analyzed on a personal computer as well as printed out as text data from a personal computer.



NIRO is a registered trademark of Hamamatsu Photonics K.K.

Product and software package names noted in this documentation are trademarks or registered trademarks of their respective manufacturers.

Subject to local technical requirements and regulations, availability of products included in this promotional material may vary. Please consult your local sales representative.

Information furnished by HAMAMATSU is believed to be reliable. However, no responsibility is assumed for possible inaccuracies or omissions.

Specifications and external appearance are subject to change without notice.

© 2010 Hamamatsu Photonics K.K.



HAMAMATSU

Homepage Address <http://www.hamamatsu.com>

HAMAMATSU PHOTONICS K.K., Systems Division
812 Joko-cho, Higashi-ku, Hamamatsu City, 431-3196, Japan, Telephone: (81)53-431-0124, Fax: (81)53-433-8031, E-mail: export@sys.hpk.co.jp

U.S.A. and Canada: Hamamatsu Corporation, 360 Foothill Road, Bridgewater, N.J. 08807-0910, U.S.A., Telephone: (1) 908-231-0960, Fax: (1) 908-231-0852, E-mail: usa@hamamatsu.com

Germany: Hamamatsu Photonics Deutschland GmbH, Arzbergerstr. 10, D-82211 Herrsching am Ammersee, Germany, Telephone: (49)8152-375-0, Fax: (49)8152-2658, E-mail: info@hamamatsu.de

France: Hamamatsu Photonics France S.A.R.L.: 19, Rue du Saule Trapu, Parc du Moulin de Massy, 91882 Massy Cedex, France, Telephone: (33)1 69 53 71 00, Fax: (33)1 69 53 71 10, E-mail: infos@hamamatsu.fr

United Kingdom: Hamamatsu Photonics UK Limited: 2 Howard Court, 10 Tewin Road, Welwyn Garden City, Hertfordshire, AL7 1BW, U.K., Telephone: (44) 1707-294888, Fax: (44) 1707-325777, E-mail: info@hamamatsu.co.uk

North Europe: Hamamatsu Photonics Norden AB: Smidesvägen 12, SE-171-41 Solna, Sweden, Telephone: (46)8-509-031-00, Fax: (46)8-509-031-01, E-mail: info@hamamatsu.se

Italy: Hamamatsu Photonics Italia S.R.L.: Strada della Moia, 1 int.6-20020 Arese (Milano), Italy, Telephone: (39)02-935 81 733, Fax: (39)02-935 81 741, E-mail: info@hamamatsu.it

Cat.No.SMES0014E01
JUL/2010 HPK
Created in Japan